

Stylish on the Outside. Smart on the Inside.



THE F&P ICON™ PREMO DELIVERS:

- ThermoSmart™ Technology for more humidity and comfort
- InfoSmart™ Technologies for flexible data communication options

PLUS

- Contemporary and stylish design for bedroom appeal
- Forward facing clock with customizable AlarmTunes™
- Compact footprint including humidifier and power supply
- Combined smart technologies.



Stylish on the Outside. Smart on the Inside.


The F&P ICON™ has been designed from the outside-in to answer the strong call for a contemporary, compact CPAP that blends seamlessly into any bedroom environment, helping to break down some of the traditional resistance associated with CPAP therapy. Inside the F&P ICON™ is a comprehensive set of intelligent technologies, integrated into one stylish solution to deliver new levels of therapy, while effectively meeting your business needs.

ThermoSmart™

More humidity. More comfort.



ThermoSmart™ Technology includes an advanced humidification algorithm and a heated breathing tube, which both work to provide optimal humidity and reduce condensation. Patients will benefit from higher levels of humidity¹ for improved comfort^{2,3} and a clinically proven better night's sleep.³

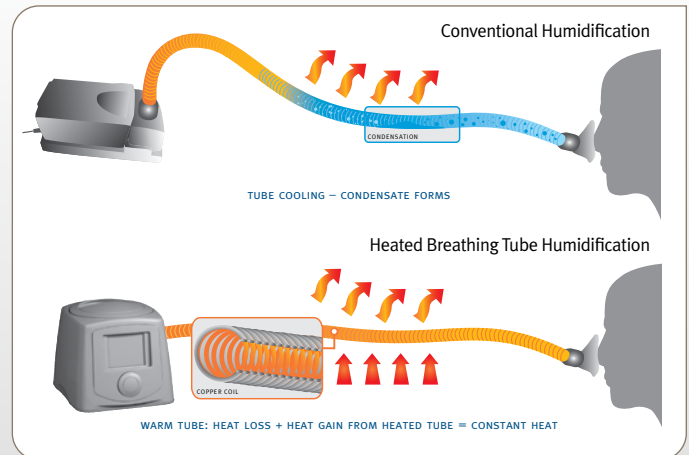

InfoSmart™ Technologies

More communication. More results.



Retrieving compliance data in a timely and efficient manner is crucial to improving patient outcomes and achieving business success. The F&P ICON™ provides a range of solutions to suit varied business needs.

Data Communication

- InfoSend™ View: The patient can relay sleep data from the display screen of the device over the telephone, which can then be verified using the F&P InfoSmart™ software
- SmartStick™ (Net): The functionality of the SmartStick™ and SmartStick™ Net have been combined into one, so the choice of data retrieval becomes even more simple and cost effective. The physical SmartStick™ can be sent via the mail for data download, or the sleep data can be sent electronically via the web.

Data Reporting

Treatment compliance and efficacy reporting is a key factor for comprehensive patient follow-up and thus therapy outcomes. F&P InfoSmart™ Software produces comprehensive charting and reporting of patient data including:

- A choice of 30, 60 or 90 days of summary data
- 7 days of detailed data
- Cumulative summary data from first use



FOCUS ON FEATURES



Elbow

Used to connect a ThermoSmart™ or standard Breathing Tube to the device. The Elbow also rotates for a preferred position.

Water Chamber

One-piece chamber. Holds 420 ml of water and is dishwasher.



Chamber Lid

Must be correctly fitted for optimal pressure delivery and can also be used as a handle to carry the device.

Menu System

The Home Screen with a clock is the default display. Each point of the clock-face represents the available settings.



SmartStick™ / USB Port

The USB Port houses the SmartStick™ to record patient compliance, efficacy and flow data.

SmartDial™

Intuitive controls allow one-touch navigation of the Menu System, using a Press, Turn, Press principle.



Air Filter

Ensures the air is filtered from particles and dust.



PREMO CPAP | FEATURES

THE F&P ICON™ PREMO

| INCLUDED WITH THE F&P ICON™ | REPLACEMENT PARTS | PRODUCT SPECIFICATIONS | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------|------------|-------------|------------|-----------------|------------|-------|------------|-------------|------------|-----------------------------|------------|----------------|----------|------------|------------|--------------|------------|------------------|------------|---------------|------------|-----------------------|------------|-----------------|------------|---------------------|--|------------|-----------|-------------|-----------|------------------------|------------|-------------------|------------|-----------------------|------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------|---------------------|------------------------|-----------|---------------------|-----------------|-----------|---------------------|----------|
| <ul style="list-style-type: none"> Premo CPAP 1 x Carry-case* 1 x Water Chamber 1 x Chamber Lid 1 x Breathing Tube** 1 x SmartStick™ 1 x SmartStick™ Mailer 1 x Use and Care Manual 1 x spare Air Filter | <table border="0"> <tr> <td>Water Chamber</td> <td>900ICON200</td> </tr> <tr> <td>SmartStick™</td> <td>900ICON202</td> </tr> <tr> <td>SmartStick™ Net</td> <td>900ICON203</td> </tr> <tr> <td>Elbow</td> <td>900ICON204</td> </tr> <tr> <td>Outlet Seal</td> <td>900ICON206</td> </tr> <tr> <td>ThermoSmart™ Breathing Tube</td> <td>900ICON208</td> </tr> <tr> <td>Breathing Tube</td> <td>900HC221</td> </tr> <tr> <td>Air Filter</td> <td>900ICON503</td> </tr> <tr> <td>Filter Grill</td> <td>900ICON218</td> </tr> <tr> <td>Silver Gloss Lid</td> <td>900ICON214</td> </tr> <tr> <td>Matt Gray Lid</td> <td>900ICON216</td> </tr> <tr> <td>Protective Carry-case</td> <td>900ICON308</td> </tr> <tr> <td>Soft Carry-case</td> <td>900ICON315</td> </tr> <tr> <td>Use and Care Manual</td> <td></td> </tr> <tr> <td> 3-Language</td> <td>185046221</td> </tr> <tr> <td> 17-Language</td> <td>185045354</td> </tr> <tr> <td>Healthcare Provider CD</td> <td>900ICON110</td> </tr> <tr> <td>F&P InfoSmart™ CD</td> <td>900ICON100</td> </tr> <tr> <td>SmartStick™ Studio CD</td> <td>900ICON112</td> </tr> </table> | Water Chamber | 900ICON200 | SmartStick™ | 900ICON202 | SmartStick™ Net | 900ICON203 | Elbow | 900ICON204 | Outlet Seal | 900ICON206 | ThermoSmart™ Breathing Tube | 900ICON208 | Breathing Tube | 900HC221 | Air Filter | 900ICON503 | Filter Grill | 900ICON218 | Silver Gloss Lid | 900ICON214 | Matt Gray Lid | 900ICON216 | Protective Carry-case | 900ICON308 | Soft Carry-case | 900ICON315 | Use and Care Manual | | 3-Language | 185046221 | 17-Language | 185045354 | Healthcare Provider CD | 900ICON110 | F&P InfoSmart™ CD | 900ICON100 | SmartStick™ Studio CD | 900ICON112 | <p>Dimensions 6.3 x 6.7 x 8.7 in (160 x 170 x 220 mm)</p> <p>Weight 4.8 lb (2.2 kg) Packed Weight: 8.7 lb (4.0 kg)</p> <p>Performance Pressure Range: 4 to 20 cmH₂O</p> <p>Noise Level < 30 dBA</p> <p>Standards Compliance IEC60601-1: 1998, AS/NZS 3200.1.0: 1998, EN60601-1: 1990 + A1 - A2</p> <p>SmartStick™ Data Recording The SmartStick™ records up to 365 days of summary compliance and efficacy data, 30 days of detailed efficacy data (AHI, Pressure, Leak), as well as cumulative data from first use.</p> <p>Electrical Ratings</p> <table border="1"> <thead> <tr> <th>RATED SUPPLY VOLTAGE</th> <th>RATED CURRENT INPUT</th> <th>RATED SUPPLY FREQUENCY</th> </tr> </thead> <tbody> <tr> <td>100-115 V</td> <td>1.27 A (1.43 A Max)</td> <td>50-60 Hz/400 Hz</td> </tr> <tr> <td>220-240 V</td> <td>1.07 A (1.21 A Max)</td> <td>50-60 Hz</td> </tr> </tbody> </table> | RATED SUPPLY VOLTAGE | RATED CURRENT INPUT | RATED SUPPLY FREQUENCY | 100-115 V | 1.27 A (1.43 A Max) | 50-60 Hz/400 Hz | 220-240 V | 1.07 A (1.21 A Max) | 50-60 Hz |
| Water Chamber | 900ICON200 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| SmartStick™ | 900ICON202 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| SmartStick™ Net | 900ICON203 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Elbow | 900ICON204 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Outlet Seal | 900ICON206 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| ThermoSmart™ Breathing Tube | 900ICON208 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Breathing Tube | 900HC221 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Air Filter | 900ICON503 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Filter Grill | 900ICON218 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Silver Gloss Lid | 900ICON214 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Matt Gray Lid | 900ICON216 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Protective Carry-case | 900ICON308 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Soft Carry-case | 900ICON315 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Use and Care Manual | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 3-Language | 185046221 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 17-Language | 185045354 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Healthcare Provider CD | 900ICON110 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| F&P InfoSmart™ CD | 900ICON100 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| SmartStick™ Studio CD | 900ICON112 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| RATED SUPPLY VOLTAGE | RATED CURRENT INPUT | RATED SUPPLY FREQUENCY | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 100-115 V | 1.27 A (1.43 A Max) | 50-60 Hz/400 Hz | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| 220-240 V | 1.07 A (1.21 A Max) | 50-60 Hz | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |

F&P ICON™ Models and Features Matrix

| PERFORMANCE FEATURES | Auto | | Premo | Novo |
|---------------------------|-----------|-----------|-----------|---------|
| | ICONAAX | ICONABX | ICONPBX | ICONNAX |
| Fully integrated | • | • | • | • |
| ThermoSmart™ Technology** | • | • | • | • |
| Auto-Adjusting Pressure | • | • | | |
| Efficacy Reporting | • | • | • | |
| Compliance Reporting | • | • | • | • |
| SmartStick™ | • | • | • | • |
| SensAwake™ | • | | | |
| Proportional Ramp | • | • | • | • |
| Auto-Altitude Adjusting | Automatic | Automatic | Automatic | Manual |
| Leak Compensation | • | • | • | |
| Clock and AlarmTunes™ | • | • | • | • |
| InfoSmart™ Technologies | Advanced | Advanced | Advanced | Basic |

NOTE: The X denotes country variant

* Part may vary depending on country

** In some countries the ThermoSmart™ Breathing Tube needs to be purchased as an accessory to activate ThermoSmart™ Technology

www.fphcare.com

Manufacturer Fisher & Paykel Healthcare Ltd, 15 Maurice Paykel Place, East Tamaki, Auckland 2013, PO Box 14 348, Panmure, Auckland 1741, New Zealand Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.co.nz Web: www.fphcare.com **Australia** Tel: +61 3 9879 5022 Fax: +61 3 9879 5232 **Austria** Tel: 0800 29 31 23 Fax: 0800 29 31 22 **Benelux** Tel: +31 40 216 3555 Fax: +31 40 216 3554 **China** Tel: +86 20 3205 3486 Fax: +86 20 3205 2132 **France** Tel: +33 1 6446 5201 Fax: +33 1 6446 5221 **Germany** Tel: +49 7181 98599 0 Fax: +49 7181 98599 66 **India** Tel: +91 80 4284 4000 Fax: +91 80 4123 6044 **Irish Republic** Tel: 1800 409 011 **Italy** Tel: +39 06 7839 2939 Fax: +39 06 7814 7709 **Japan** Tel: +81 3 3661 7205 Fax: +81 3 3661 7206 **Northern Ireland** Tel: +0800 132 189 **Spain** Tel: +34 902 013 346 Fax: +34 902 013 379 **Sweden** Tel: +46 8 564 76 680 Fax: +46 8 36 63 10 **Switzerland** Tel: 0800 83 47 63 Fax: 0800 83 47 54 **Taiwan** Tel: +886 2 8751 1739 Fax: +886 2 8751 5625 **Turkey** Tel: +90 312 354 34 12 Fax: +90 312 354 31 01 **UK (EU Authorised Representative)** Fisher & Paykel Healthcare Ltd, Unit 16, Cordwallis Park, Clivemont Road, Maidenhead, Berkshire SL6 7BU, UK Tel: +44 1628 626 136 Fax: +44 1628 626 146 **USA** Tel: +1 800 446 3908 or +1 949 453 4000 Fax: +1 949 453 4001